




# Long term visual results of Prima chip in patients with geographic atrophy



R. Hornig, Clinical Affairs Director, Pixium Vision

- Forward Looking Statements** This document contains information on Pixium Vision's markets and competitive position, and more specifically, on the size of its potential markets. This information has been drawn from various sources or from the companies own estimates. Investors should not base their investment decision on this information. This document also contains certain forward-looking statements. These statements are not a guarantee of the Company's future performance. These forward-looking statements relate to the Company's future prospects, developments and marketing strategy and are based on analyses of earnings forecasts and estimates of amounts not yet determinable. Forward-looking statements are subject to a variety of risks and uncertainties as they relate to future events and are dependent on circumstances that may or may not materialize in the future. Pixium Vision draws your attention to the fact that as forward-looking statements cannot under any circumstance be construed as a guarantee of the Company's future performance and that the Company's actual financial position, results and cash flow, as well as the trends in the sector in which the Company operates were to conform to the forward-looking statements contained in this document. Furthermore, even if Pixium Vision's financial position, results, cash flow and developments in the sector in which the Company operates were to conform to the forward-looking statements contained in this document, such results or developments cannot be construed as a reliable indication of the Company's future results or developments. The Company does not undertake any obligation to update or to confirm projections or estimates made by analysts or to make public any correction to any prospective information in order to reflect an event or circumstance that may occur after the date of this presentation. A description of those events that may have a material adverse effect on the business, financial position or results of Pixium Vision, or on its ability to meet its targets, appears in the sections "Risk Factors" of its "Document de Base" filed with the French Autorité des Marchés Financiers. By attending this presentation or accepting this document, you agree to be bound by the foregoing restrictions set out above.
- Purpose:** To evaluate the long term central prosthetic vision with the photovoltaic subretinal implant activated by augmented-reality glasses and simultaneous perception of the natural peripheral vision in patients with geographic atrophy.
- Methods:** Five patients with visual acuity  $\leq 20/400$  due to geographic atrophy of at least 3 optic disc diameters and no foveal vision have been implanted with a wireless photovoltaic chip (PRIMA, Pixium Vision) of 2x2mm in size, 30 $\mu$ m in thickness, containing 378 pixels of 100 $\mu$ m in width. Each pixel in the implant converts pulsed near-infrared light (880nm) projected from video glasses into electric current to stimulate the nearby neurons in the inner nuclear layer of the retina. Prosthetic acuity was assessed using electronic magnification. Simultaneous perception of central prosthetic and peripheral natural vision was evaluated under room lighting.
- Results:** In all patients, the chip implanted under the macula remained stable and functional. One patient passed away after 18 months, and the 4 others have a follow-up extending now 3 years. No decrease in natural eccentric visual acuity was observed in any of the study eyes. All 5 patients could perceive white-yellow patterns with adjustable brightness, in retinotopically correct locations within previous scotoma during the follow-up period. With electronic magnification visual acuity of 20/63 was reached. Under room lighting, patients could simultaneously use prosthetic central vision and the remaining peripheral vision in the operated eye and in the fellow eye.
- Conclusions:** Wireless chip PRIMA implanted under the atrophic macula in patients with geographic atrophy remains stable and functional during the 3 years of follow-up. The implant provides central visual perception with acuity close to the single pixel size of the photovoltaic array. Augmented reality glasses enable simultaneous perception of the central prosthetic and natural peripheral vision under room lighting, while electronic zoom provides significantly higher resolution.

## PRIMA mechanism of action

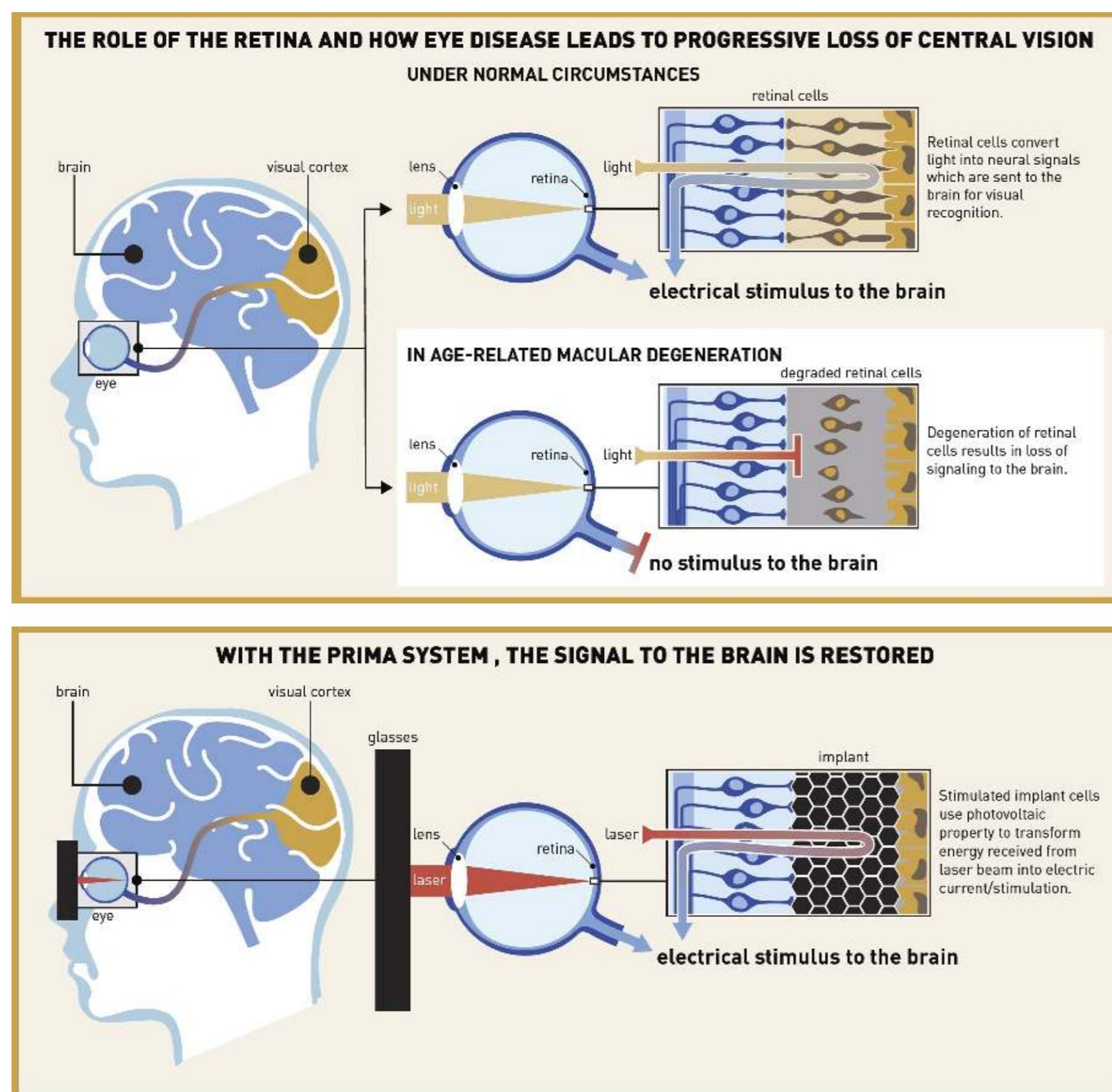

## The PRIMA retinal implant system



Illustrations for information purposes – not indicative of actual size or clinical outcome

## Implant and Delivery System Design



### Implant features:

- Implant size: 2 mm x 2 mm x 30  $\mu$ m
- 378 pixels
- Pixel size: 100  $\mu$ m

## PRIMA, Human Feasibility Studies in patients suffering from Advanced dry-AMD

**FEASIBILITY STUDIES**

**Feasibility Study of Compensation for Blindness With the PRIMA System in Patients With Dry Age Related Macular Degeneration**

**CRITERIA**

- 60 years or older
- Visual Acuity of LogMAR 1.3 (20/400) or worse
- No foveal perception

**ENDPOINTS**

Visual perception with PRIMA implant, visual acuity and safety assessment after implantation.

**PRIMA-FS** ClinicalTrials.gov  
5 patients implanted in France  
Recruitment closed  
Positive interim data available

**PRIMA-FS-US** ClinicalTrials.gov  
5 patients in USA – Recruitment open  
Safety assessment at 12 month after implantation

**CLINICAL SITES**

Dr. Y. Le Mer  
Dr. J-F Girmens  
Fondation Ophtalmologique Adolphe de Rothschild

Pr. J. Sahel  
Dr. J. Martel  
UPMC

Pr. B. Lam,  
Pr. N. Gregori  
Pr J.-M. Parel  
Bascom Palmer Eye Institute  
UNIVERSITY OF MIAMI HEALTH SYSTEM

## Feasibility Study in France


### Feasibility Study of Compensation for Blindness With the PRIMA System in Patients With Dry Age Related Macular Degeneration

The primary objective of the feasibility trial is to demonstrate that the subjects have visual perception elicited by the PRIMA System.

### The key-inclusion criteria imply that the patient:

- Is **60 years** or older at the date of inclusion;
- Has a confirmed diagnosis of **advanced dry age related macular degeneration** with an atrophy size of at least 3 optic disc diameters;
- Has best corrected visual acuity in the study eye of **logMAR 1.3 (20/400) or worse** measured by ETDRS;
- Has **no foveal perception** measured by micro-perimetry in the study eye;
- Has a **refraction of study eye between -3 and +4** (limits included) for patients with IOL (there is no refraction criteria for patients with natural lens);
- Understands and accepts the obligation to present for all schedule follow-up visits.

## PRIMA Feasibility Anatomical Results

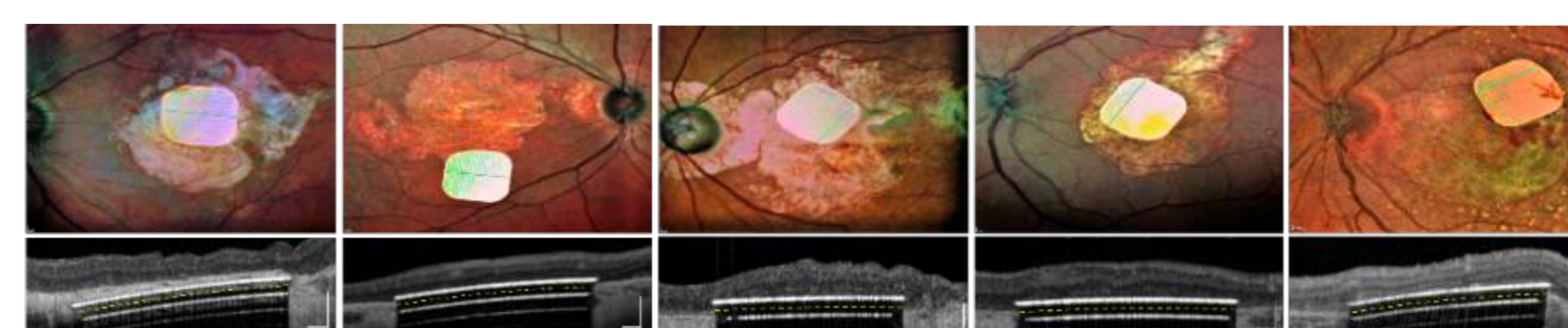


\*The implant appears in OCT twice thicker due to higher refractive index of Si compared to that of the retinal tissue. Dashed line indicates position of the back side of the implant resting on Bruch's membrane.

Implant location in the macula	Intra-chorooidal	Central subretinal	Central subretinal	Off-center subretinal	Central subretinal
Average distance from implant to bottom of INL at 3, 6, 12 months ( $\mu$ m)	127	34	43	37	51
	103	28	35	38	39
	138	39	39	35	37


Palanker et al., Photovoltaic Restoration of Central Vision in Atrophic Age-Related Macular Degeneration, Ophthalmology 2020

## Natural visual acuity (ETDRS) of study eye in Feasibility study

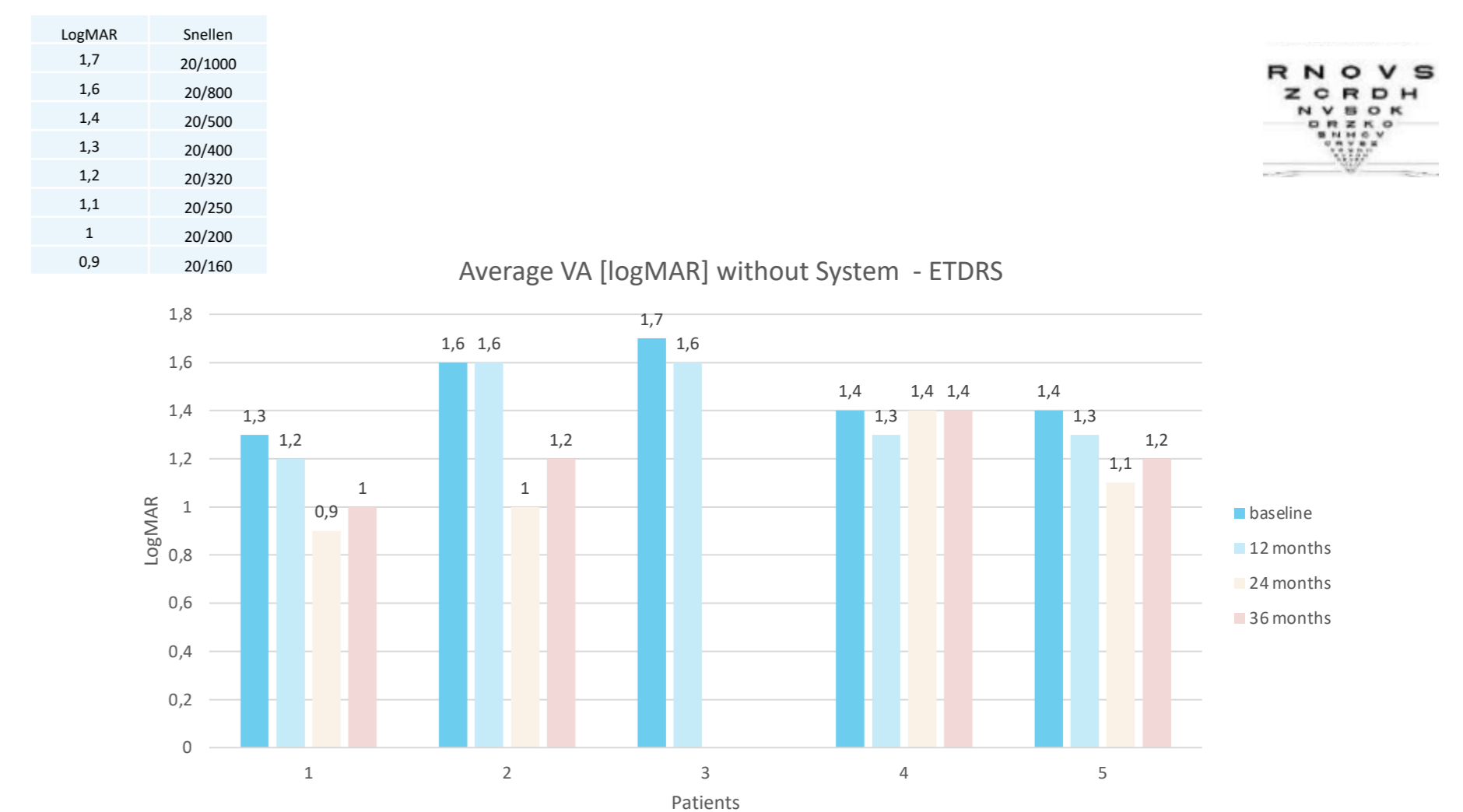

PRIMA could be safely implanted under the atrophic macula while preserving the residual natural peripheral visual acuity, measured via ETDRS in the Feasibility study

## PRIMA Feasibility visual acuity tests

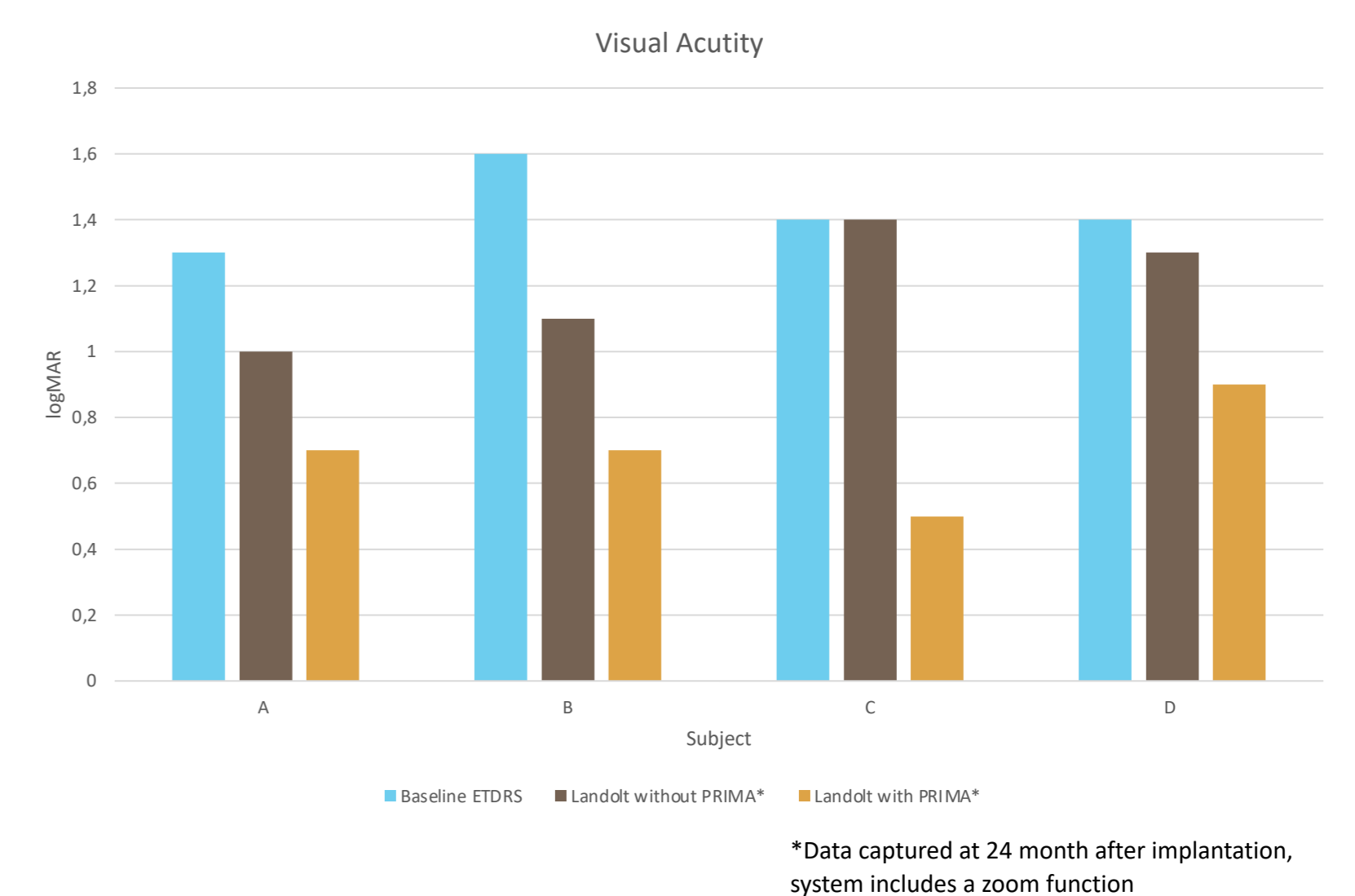

## PRIMA Feasibility daily living observations



Single case reports may not be representative for other patients using the PRIMA System

## Feasibility studies : Conclusions

- Subretinal implantation of the wireless PRIMA microchip in atrophic dry AMD is feasible.
- No decrease in residual natural vision compared to pre-operative visual acuity.
- ALL 5 patients reached the functional primary end point of the feasibility study (visual light perception in the former scotoma).
- PRIMA provided the best prosthetic visual acuity to date
- Most subjects are able to read single letters and small words
- Some patients are able to read words and small sentences in different light conditions (Outdoors and ambient light)